


SHI-174 : SIEVE SHAKER (ROTOP MOTORIZED)

Our RoTap Sieve Shaker machines have been engineered with reliable, rugged steel and alloy materials. The mechanism is fully immersed in oil tank. These units are fitted with heavy duty **electric motor (1/4 HP single phase 220-230 V 50 Hz.)** that enables the machine to provide efficient circular and tapping motion given by hammering at the top

Our standard model is designed to accept up to 8" diameter sieves and also accommodates

Regardless of diameter, any one material sample can be tested on a series of 7 full height (2") sieves at one time, plus bottom pan and top cover lid.

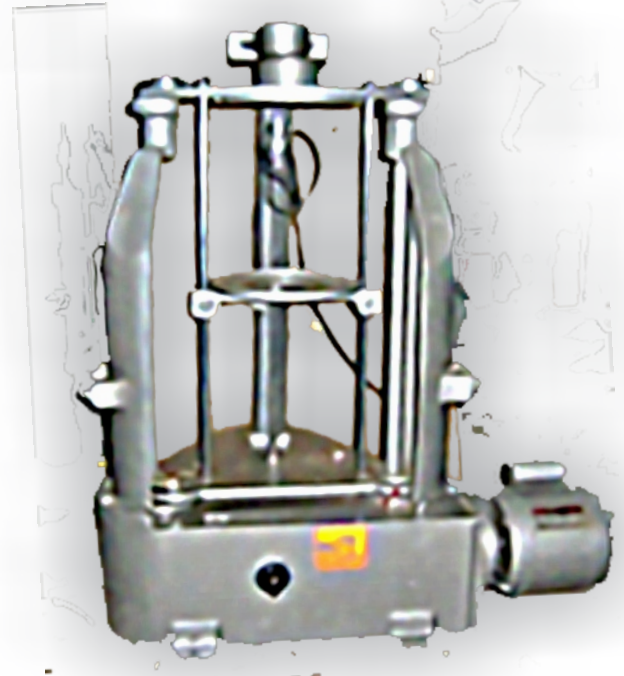

Model	SHI-174
Capacity	7sieves of 200mm dia and 50mm height with top lid &pan. 7sieves of 300mm dia and 50mm height with top lid &pan.(Optional)
Construction	Cast Iron, Painted Mild
Circular motion	250RPM
Tapping	150 minimum
Motor	¼ HP Geared mechanism
Timer	0 - 60 Minutes
Timer Type	Analog/digital (Optional)
Standard	ASTM
Power Supply	220-240 Volts 50Hz.
Certifications	CE and ISO

SHIVAM INSTRUMENTS
DELHI-110052
9312838091, 9811362676



mverma746@gmail.com | info@shivamlabinstruments.com
www.shivamlabinstruments.com